

# Blood Typing Manager

All your results, secured at your fingertips.  
Easy. Efficient. Secure.

An innovative compliant data management tool to centralize multi-instrument and multi-site data handling securely from a single interface, facilitating implementation, thus reducing IT efforts and budget.



## Easy implementation

- Benefit from the **flexible and very easy installation**, which can be performed in a physical server or virtual environment
- Optimize application set-up for one or multiple laboratories from a **single interface**, reducing resources needed to manage and maintain the system
- Achieve multi-site and multi-instrument connectivity with this **web-based software** compatible with Windows 10 and requiring no special client workstations



## Smooth connectivity

- Efficiently connect your IH instruments to LIS with the assistance of **Grifols subject matter experts**
- Achieve a **smooth and unique** bidirectional connection to LIS
- Improve work efficiency of your lab or network with the ability to interface to **multiple LIS**




## Secured data

- Have confidence knowing that the confidentiality, integrity, and availability of your data are in good hands. Grifols designs products with **cybersecurity** in mind
- Patient data is securely managed in Blood Typing Manager, complying with **data protection** and **cybersecurity standards**
- Keep your lab running smoothly with Secure Remote Connectivity, a tool for the technical support team to **remotely access** the middleware to reduce incident handling time while optimizing system uptime

TYPING

# Blood Typing Manager

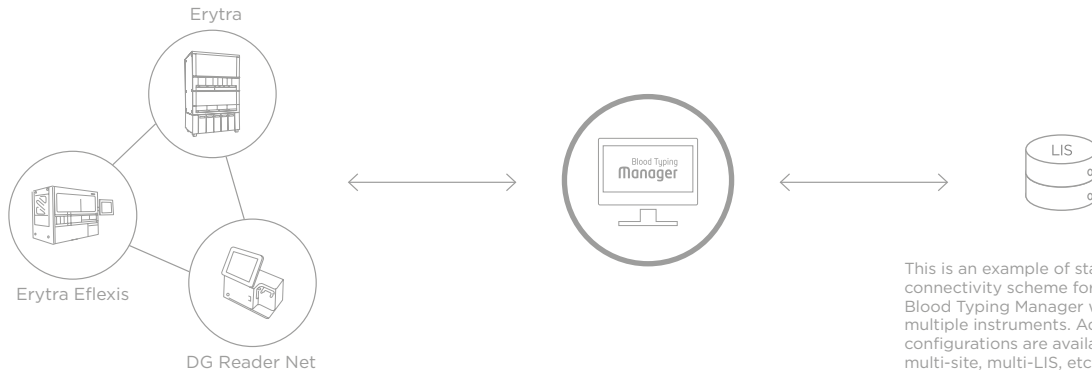

This is an example of standard connectivity scheme for Blood Typing Manager with multiple instruments. Additional configurations are available: multi-site, multi-LIS, etc.

## Technical specifications

### INSTALLATION

Component	Server Computer Specification	Client Computer Specification
<b>CPU</b>	Intel 4 cores 3GHz or faster	Intel Celeron or faster
<b>Memory</b>	8 GB or more	4 GB or more
<b>Hard drive</b>	250 GB or more, 7200 rpm or faster	250 GB or more
<b>Operating system</b>	Windows 10 v1607 (or higher), or Windows Server 2016 (or higher)*	Not applicable
<b>Web browser</b>	Not applicable	Google Chrome and Microsoft Edge

### ORDERING INFORMATION

Code	Description
<b>21040401</b>	Blood Typing Manager Core Module
<b>21040402</b>	Blood Typing Manager Multi-site Plus (> 5 sites)
<b>21040403</b>	Blood Typing Manager Antibody Identification Module

\*Windows 10 is recommended for environments of up to 4 concurrent users and 5 connected instruments. In other environments, Windows Server 2016 or above is recommended.

### SOFTWARE CONNECTIVITY

<b>LIS connectivity</b>	Bidirectional connectivity to one or more LIS
<b>LIS protocol</b>	Serial port: HL7, File sharing: HL7, TCP/IP: HL7
<b>Instrument connectivity</b>	File sharing or TCP/IP connectivity to unlimited number of IH Grifols instruments
<b>IH Grifols instruments compatibility</b>	Erytra, Erytra Eflexis, DG Reader Net
<b>Sites connectivity</b>	Multiple laboratory connectivity sharing a unique database
<b>Client PC connectivity</b>	Web access (HTTPS)
<b>Technical support and remote connectivity</b>	Compatible with Grifols Secure Remote Support tool that provides remote desktop and file transfer capabilities for troubleshooting for technical support

### STANDARDS

<b>Regulatory standards</b>	CE Mark, FDA 510k; ISO 13485 & ISO 14971; IEC 62304 & IEC 62366-1, UKCA mark according to UK Medical Devices Regulation 2002
<b>Cybersecurity standards</b>	AAMI TIR57: Principles for medical device security
<b>Data protection standards</b>	HIPAA (Health Insurance Portability and Accountability Act), 45 CFR Parts 160, 164 Regulation (EU) 2016/679 (General Data Protection Regulation)

DIAGNOSTIC  
**Support & Services**

### With you every step of the way

We deliver solutions to facilitate high-quality patient care with a dedicated and responsive team of professionals.

**GRIFOLS**

**Diagnostic Grifols, S.A.**  
Pg. Fluvial, 24  
08015 Parets del Vallès  
Barcelona, Spain  
Tel: +34 935 710 400

**Grifols International, S.A.**  
Parc Empresarial Can Sant Joan  
Av. de la Generalitat, 152-158  
08174 Sant Cugat del Vallès  
Barcelona, Spain  
Tel: +34 935 710 500